

Software release notes

1. AFFECTED VERSION

NOVAscope software version **1.15P** is intended for use on **Windows 7** systems.



In case the NOVA runs on SW 1.9 or older, it is recommended to upgrade to software 1.10C first and followingly proceed with upgrade 1.15P.



NOVAscope software version 1.15P is **not** compatible with Windows 10 systems.




The PC application to review data is the same for Windows 7 and Windows 10 systems.

For further details, please contact support@finapres.com

2. HOW TO CHECK THE WINDOWS VERSION ON YOUR NOVA

Windows version can be checked on the NOVA following the described instructions:

Main Menu > System Settings > Software > Current Version

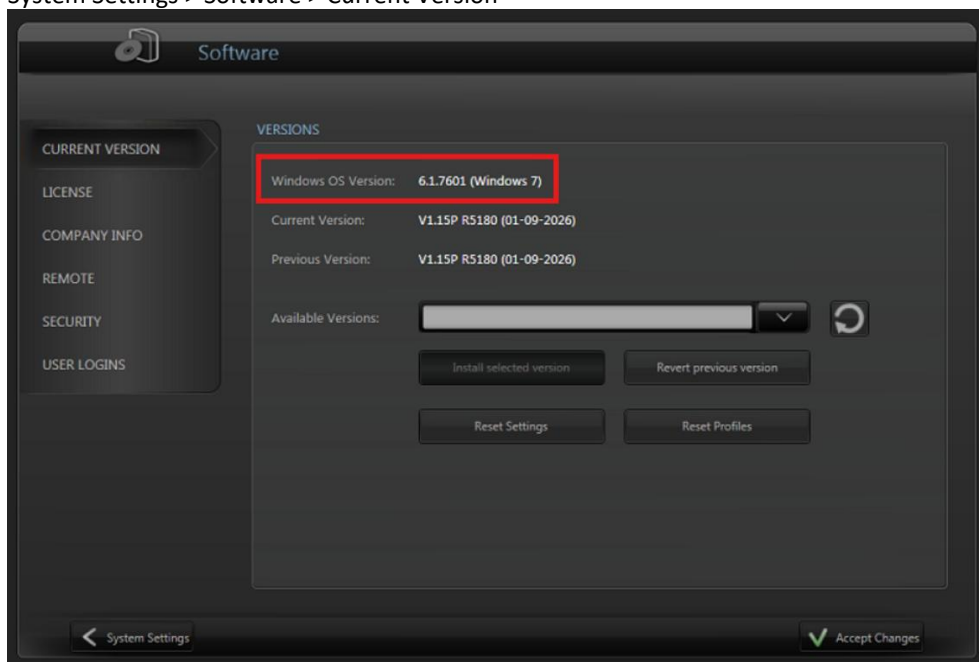

3. NEW FEATURES

No new features

4. IMPROVEMENTS

- **Improved visibility of selected button:** When a button is selected (e.g. zoom in button during measurement review), it is now greyed and more visible by the user.

5. BUG FIXES

- **Arm cuff values:** When an arm cuff measurement is performed during the measurement, the numerical values will not disappear after 1 minute.
- **iPad app - Valsalva feedback bar:** The feedback bar shown on the iPad during Valsalva maneuver has been restored.
- **Encryption tab:** Encryption tab is now accessible from Software System Settings (the option is not grayed out anymore).
- **Keyboard special character:** When special character combinations are typed, nothing happens.
- **Prepare for measurement window:**
 - Removed the dot option in the height field;
 - Enabled “Next” and “Back” button when age field is selected.